

PowderGuard+™

Ductless Fume Cupboard



The PowderGuard+™ ductless fume cupboard provides a safe, contained environment for powder weighing and other compounding processes in conjunction with noxious chemicals. Examination and manipulation of toxic compounds can be achieved whilst protecting the operator, using a state-of-the-art HEPA filtration combination. Australian designed and manufactured, the PowderGuard+ design promotes durability and safety, and requires no ductwork and minimal installation.

Filter replacement and airflow are carefully monitored through an operator control system, which will alert the user of changes when necessary. HEPA and pre-filter replacement can be undertaken safely and confidently by the operator, and our technicians can provide further assistance and servicing. To find out more about the PowderGuard+™, speak to one of our specialists today about your specific applications.




Features

- Ideal for compounding of toxic particulates
- Filtered, ductless design
- Designed & made in Australia
- Compact & modern for smaller environments
- Minimal installation
- Low airflow & filter replacement alarm
- Features HEPA filtration system
- Pre-filter/HEPA bag-out under containment
- Multiple sizes available
- Standard LED light

Rev 17/06



PowderGuard+™

Ductless Fume Cupboard

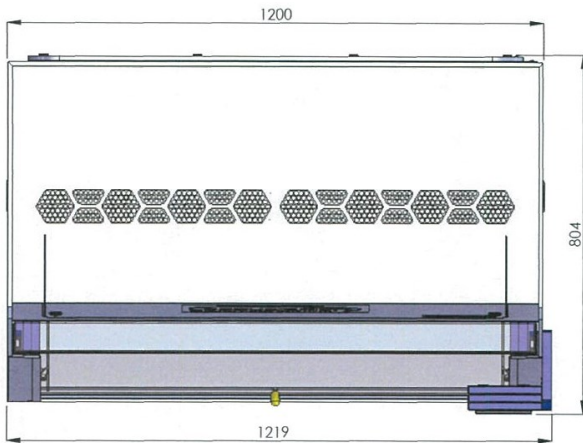

Optional Features

- Custom base inserts
- Rigid underbench frame
- Mobile under bench
- Electric height-adjustable bench

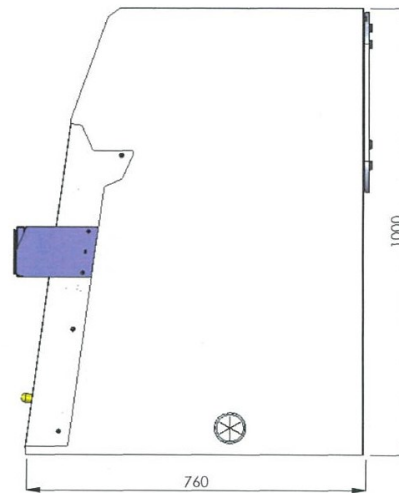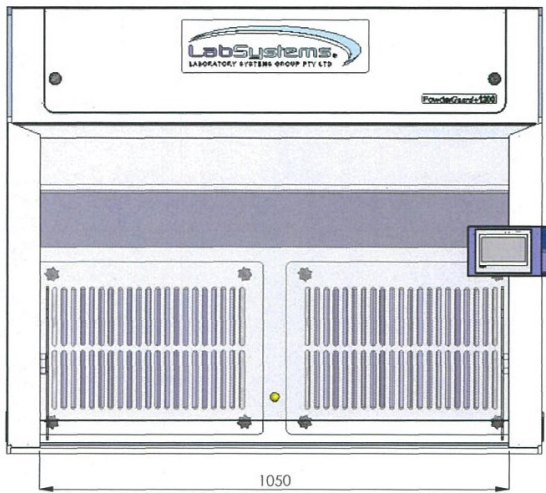

DIMENSIONS	PowderGuard + 800™	PowderGuard + 1200™
HEIGHT—EXTERNAL	1000MM	1000MM
WIDTH—EXTERNAL	819MM	1219MM
DEPTH—EXTERNAL	804MM	804MM